

For further information please contact

JohnsonDiversey
Weston Favell Centre
Northampton
NN3 8PD
Sales and Service
Tel: 01604 783505
Fax: 01604 783506

or visit our website at
www.johnsondiversev.co.uk



The JohnsonDiversey Advantage

Beyond clean, we are dedicated to being the best at simplifying the lives of our customers.

We are committed to doing this by:

- investing time to listen, understand and respond to our customers' unique cleaning and hygiene needs
- taking a personal interest to ensure the facilities they care for are consistently clean, safe and attractive
- pursuing innovation in every form likely to make the lives of our customers simpler and more profitable
- partnering with our customers to exceed their expectations every day, everywhere



ClearKlens - sterile alcohol range

the vision of JohnsonDiversey

Best solutions and service for you everyday, everywhere - defines our focus and supports our principles.



JohnsonDiversey is proud to introduce their new sterile ready to use alcohol disinfectants into the ClearKlens product range, designed specifically for the needs of your industry.

giving you assurance that:

- formulations are protected, minimising the need for revalidation
- disinfectants fully supported through the Biocides Product Directive (BPD)
- products are manufactured to GMP standards
- products are supported by full technical documentation
- disinfectants pass EN standards




clean is just the beginning

At JohnsonDiversey, clean is our promise. It is just the beginning of what we do. We go beyond clean to deliver total solutions to our customers. We are dedicated to anticipating our customers' needs and exceeding their expectations. We commit to this everyday, everywhere our customers operate. The pharmaceutical, cosmetic and toiletries (PCT) division's primary concern is the quality of your end product and we work with you to achieve it through every stage of your production process. We make sure you obtain the highest standards of hygiene for all production

equipment, safeguarding even your most sensitive production areas. We will help you do all this simply and consistently.

To enable you to maximise JohnsonDiversey's expertise in the pharmaceutical industry, we have developed a range of sterile, ready to use Alcohol disinfectant sprays for use in Pharmaceutical cleanrooms, complimenting the ClearKlens Tego 2000 and ClearKlens Tego 2001 sterile product range also available.

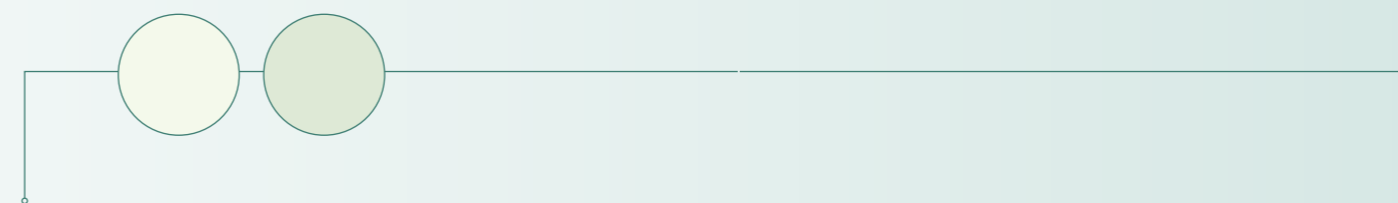

product specification

ClearKlens IPA for the pharmaceutical industry is a broad spectrum Alcohol disinfectant (70% IPA and 30% Water for Injection) in a ready to use 900ml spray format. For use on bench tops, equipment, gloved hands, instruments and utensils.



Both products are:

- Manufactured in an ISO Class 5 (Class A) environment within an ISO Class 7 (Class C) background
- 0.2 micron filtered
- Double bagged
- Gamma irradiated
- Detex irradiation indicator
- Full batch traceability
- Manufactured to a GMP standard which equates to a terminally sterilised injectable pharmaceutical
- Each batch is sterility and endotoxin tested

Documentation

In support of each batch of ClearKlens IPA and ClearKlens DE you receive:

- Certificate of Analysis
- Certificate of Irradiation
- Certificate of Sterility
- Certificate of Endotoxin Analysis

Technical Folder

A fully comprehensive Technical File for each of the Sterile Alcohol products is available. In order to make your process of validation easier, this covers:

- Full technical product specification
- Production Process
- Validation of certification
- Microbiological data

And much more...

Product	Pack Size	Code
ClearKlens IPA	6 x 900ml	7513400
ClearKlens DE	6 x 900ml	6069643

ClearKlens DE for NHS manufacturing pharmacies is a broad spectrum Alcohol disinfectant (70% Ethanol, 30% Water for Injection and 10ppm Bitrex) in a ready to use 900ml spray format. For use bench tops, equipment, gloved hands, instruments and utensils.



Packaging

ClearKlens IPA and ClearKlens DE are available in boxes of 6 x 900ml ready to use spray bottles, these contain a bag in bottle system which offers many benefits. The ready to use solution means no point of use dilution and the packaging is environmentally friendly for disposal.

Micro data

IPA and Denatured Ethanol both have a broad spectrum of kill leaving no residues and Water for Injection is used to manage the endotoxin risk. Both ClearKlens IPA and ClearKlens DE pass the following tests:

- EN 1276
- EN1650
- EN 13697

70% IPA is universally accepted by the Pharmaceutical industry and Denatured Ethanol is universally accepted by the National Health Service.

Denatured Ethanol with Bitrex - no methanol (so less toxic), no wood naphtha means controlled formulation, no variability

"Anti suck-back mechanism protects the integrity of the trigger bottle contents for months of in use shelf-life"

